

Relay 1000566

Instruction sheet

10/15 Hh



1. Safety instructions

- Do not exceed the maximal values specified for the switching voltage, switching current and the power rating specified in the technical data section.

This constitutes a fire hazard.

- Do not connect this relay to high voltage since this might break down the potential isolation.

2. Scope of delivery

- 1 Relay
- 1 Mini DIN connecting lead 8-pin, 60 cm length

3. Description

The relay is an event-triggered switch for switching on or off components that are electrically isolated from the experiment set-up. Examples of such components include light bulbs, small to medium-sized electrical motors, holding magnets and electromagnetically operated valves.

The power relay box is equipped with a set of contacts (1 N/C contact, 1 N/O contact, 1 change-over contact) with printed labels so that can be used as an N/C or an N/O switch.

When the relay coil is inactive (when no current is flowing), the N/C contact is closed.

The switch contacts consist of high-grade silver alloy and are connected to 4-mm safety sockets.

The relay is equipped with safe potential isolation between the coil and the contacts in conformance with the VDE 160 standard.

4. Technical data

Switching voltage:	250 V AC 220V DC max.
Switching current:	6A AC 0.12 A DC max.
Power rating:	1500 VA max. 1 mW min.

5. Operation

- Connect the relay to the “Digital Outputs” socket of the 3B NETlog™ equipment using the supplied mini DIN lead. It responds to digital output A.
- Preferably use 4-mm safety leads for the remaining electrical connections to the experiment set-up.

6. Sample experiment

Measuring the discharge of a capacitor with the basic experiment board

1 3B NETlog™ @ 230 V	1000540
ou	
1 3B NETlog™ @ 115 V	1000539
1 3B NETlab™	1000544
1 Relay	1000566
1 Basic experiment board @ 230 V	1000573
or	
1 Basic experiment board @ 115 V	1000572
1 Set of 15 Safety Experiment Leads, 75 cm	1002843

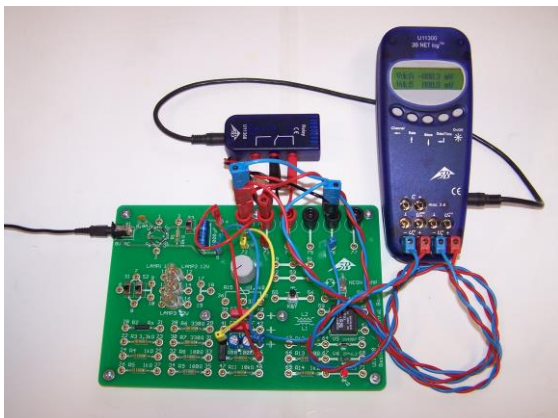


Fig. 1: Measuring the discharge of a capacitor

- Consult the instructions for the basic experiment board for the electrical connec-

tions between the relay and the experiment.

- Select one of the analog inputs (A or B) on the 3B NETlog™ module and follow the instructions for the board regarding the experiment on capacitor discharge. All necessary output value settings are specified here.
- Conduct and analyse the experiment.

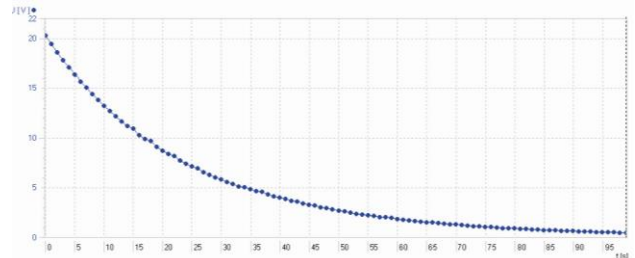


Fig. 2: Discharge curve of a capacitor