


Electroscope 1003048

Instruction Sheet

05/18 ALF

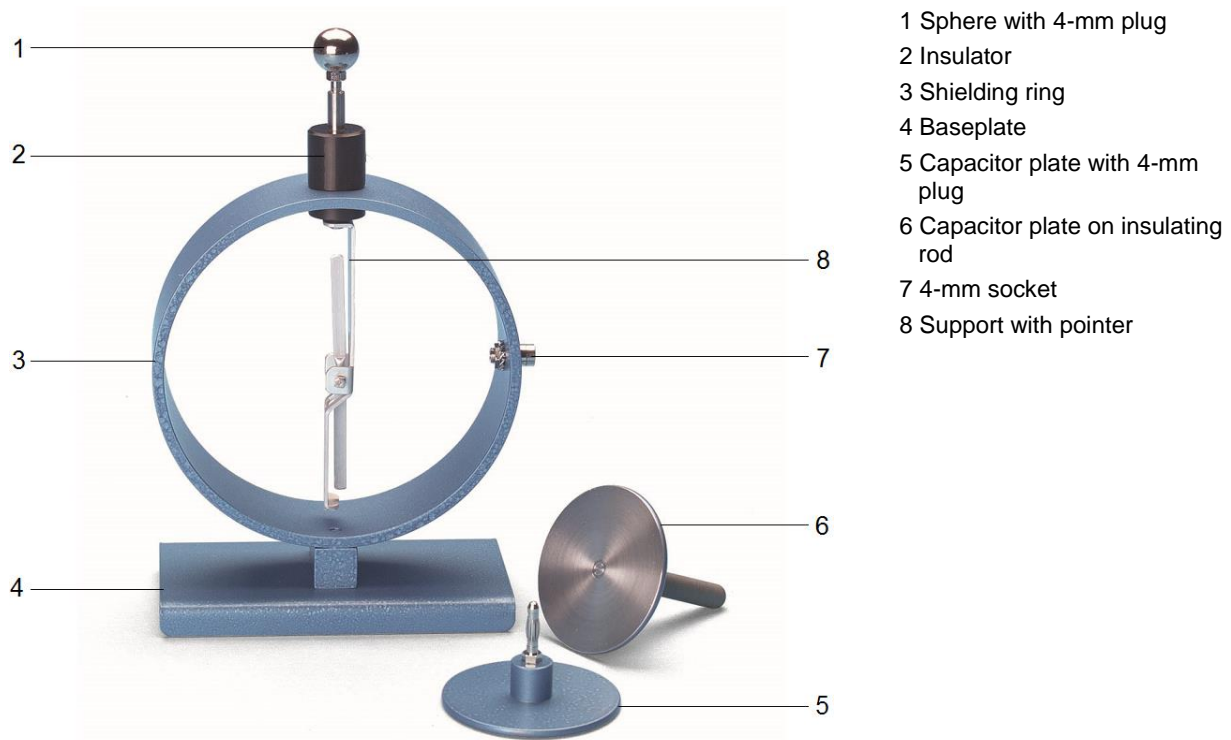

1. Description

The electroscopes is used for an introduction into the fundamentals of electrostatics and for displaying DC and AC voltages as well as demonstrating the operating principles of high voltage meters according to Braun.

It consists of a metal baseplate on which the iron shielding ring for defined field distribution is mounted. The iron ring is equipped with a 4-mm socket for earthing purposes. The support for needle with pivot bearing is attached electrically insulated to the shielding ring. The pointer is suspended asymmetrically and so its weight delivers the restoring moment. Inside the electrically conductive insulator at the upper part

of the device there is a 4-mm socket for mounting the sphere and capacitor plate.

For demonstration experiments the electroscopes is suitable for shadow projection using a point source light.

Diameter: 130 mm

2. Equipment supplied

- 1 Electroscopes
- 1 Sphere with 4-mm plug
- 1 Capacitor plate with 4-mm plug
- 1 Capacitor plate on insulating rod

3. Operation

The electroscope is charged using friction rods (e.g. 1002709) statically charged by various friction-inducing materials, or using a high-voltage power supply unit (e.g. 1003309 or 1003310).

- Make sure that the insulator is always clean and dry. If necessary use alcohol or spirits for cleaning.
- At high humidity and after transporting the unit from a cool room into a warmer one, dry the electroscope in a stream of hot air (e.g. a hair dryer).

3.1 Charging up the electroscope by touching it with a statically charged body

- Attach the sphere or capacitor plate to the electroscope.
- Rub the friction rod with the suitable material (PVC or acrylic rods e.g. with plastic foil).
- Touch the sphere or capacitor plate with the charged rod. The pointer deflects.
- Remove the friction rod, the pointer remains deflected.
- Touch the sphere with your hand. The pointer returns to normal.
- When a high-voltage power supply is used then attach the capacitor plate to the positive

pole of the power supply unit. Connect the negative pole with the earth socket on the power supply and with the earth socket on the electroscope.

- Use a charge spoon (e.g. 1002707) to draw a charge from the capacitor plate and transfer it to the electroscope

3.2 Using electrostatic induction to charge up the electroscope

- Approach but do not touch the sphere (or capacitor plate) with the statically charged friction rod. The pointer deflects.
- Remove the friction rod. The pointer returns to normal.
- Again approach the sphere with the statically charged friction rod. Once again the pointer deflects.
- Briefly touch the sphere with your finger to discharge it. The pointer deflection disappears and returns to normal.
- Now remove the friction rod. The pointer again shows deflection.

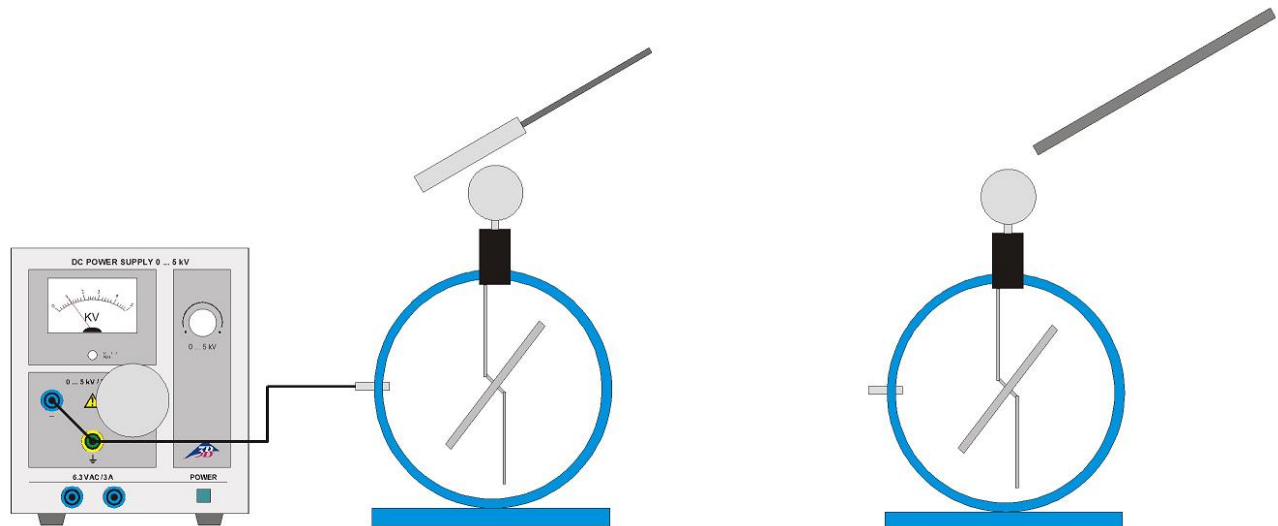


Fig. 1 Left: Charging the electroscope using the high-voltage power supply 1003309 / 1003310
Right: Charging the electroscope using electrostatic induction